

# eRx Certificate expiring 25 June 2024

The **eRx adapter** used by Bp Premier requires an update to manage the **erx.com.au certificate expiration** on **25 June 2024**.


The adapter update was included as part of the **June 2024 Bp Premier Data update**. However, some sites have reported issues with the adapter update, or have not yet run the June 2024 Data Update.

01

## Install the June 2024 Data Update

To continue sending electronic scripts via eRx Script Exchange after 25 June 2024, you must have the **June 2024 data update** installed.

To check the data update version, go to **Help > About** in Bp Premier.

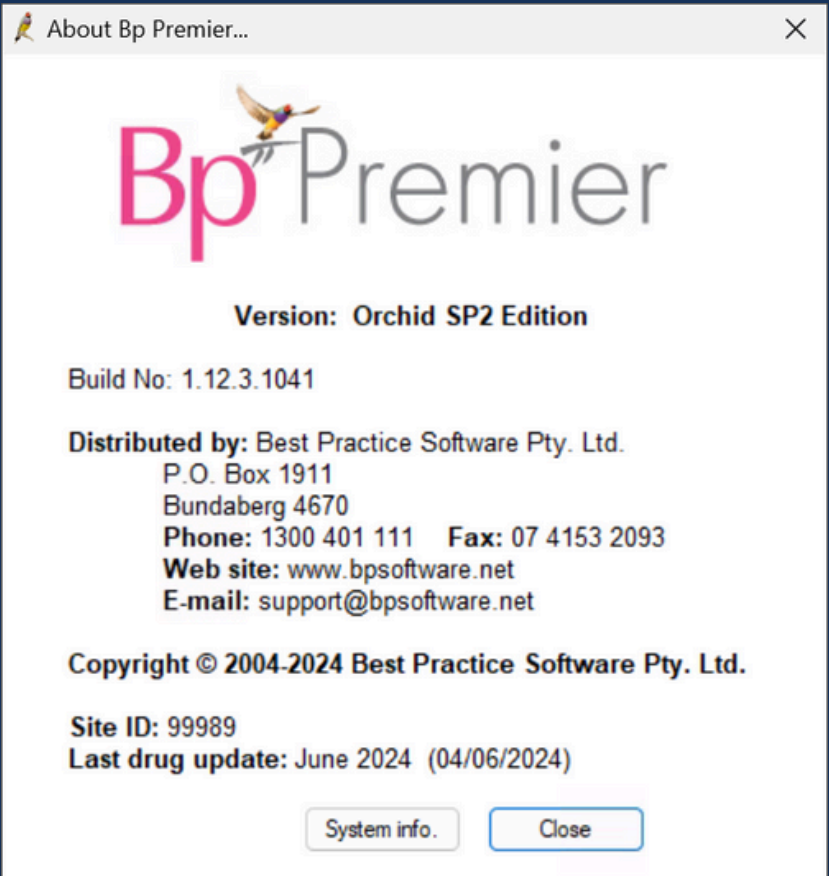

02

## Download the UpdaterWithoutDotNet\_v22.exe

Visit the [knowledge base](#) and download the [supplied .exe file](#) on the database server.

03

## Inspect the erxupdater.log in a text editor

In a Window File Explorer, browse to **C:\erx** and inspect the log file **erxupdater.log** in a text editor. Scroll to the bottom of the text file to find the **exit code**.

```
2024-06-24 14:02:42,587 [1] INFO - Services started.
2024-06-24 14:02:42,587 [1] INFO - Starting Ping test.
2024-06-24 14:02:44,165 [1] INFO - Ping test: AdapterService - 6BPX0:True
2024-06-24 14:02:44,165 [1] INFO - Ping test: Interface Service:True
2024-06-24 14:02:44,165 [1] INFO - Ping test: Business Service:True
2024-06-24 14:02:44,165 [1] INFO - Starting Ping finished.
2024-06-24 14:02:44,165 [1] INFO - Finished ...
2024-06-24 14:02:44,165 [1] INFO - Updater exiting with code 0
```

04

## Did you recieve exit code = 0

The eRx Standard Adapter has been successfully updated. No further action is required.

```
2024-06-24 14:02:42,587 [1] INFO - Services started.
2024-06-24 14:02:42,587 [1] INFO - Starting Ping test.
2024-06-24 14:02:44,165 [1] INFO - Ping test: AdapterService - 6BPX0:True
2024-06-24 14:02:44,165 [1] INFO - Ping test: Interface Service:True
2024-06-24 14:02:44,165 [1] INFO - Ping test: Business Service:True
2024-06-24 14:02:44,165 [1] INFO - Starting Ping finished.
2024-06-24 14:02:44,165 [1] INFO - Finished ...
2024-06-24 14:02:44,165 [1] INFO - Updater exiting with code 0
```

05

## Did you recieve exit code greater than 0? For example 2.

The eRx Standard Adapter has not been updated and you will need to contact Bp Premier Support.

```
--- End of inner exception stack trace ---
at ErxStandardAdapterInstaller.UpdateValidationService.PingErx()
at ErxStandardAdapterInstaller.ReleaseRunner.OnAfterExecuteRelease(IRelease release)
at ErxStandardAdapterInstaller.ReleaseRunner.ExecuteRelease(IRelease release)
at ErxStandardAdapterInstaller.ReleaseRunner.Run(String[] numbers)
at ErxStandardAdapterInstaller.Program.Main(String[] args)
2024-06-24 13:39:35,855 [1] INFO - Updater exiting with code 2
```